


Polypectomy Snares

In order to make the removal of polyps in everyday hospital routine as efficient as possible, MICRO-TECH has developed a new program of disposable polypectomy snares which perfectly combines quality with cost effectiveness.

Each oval snare is made of specially stranded wire, which convinces by its good radial force, high dimensional stability and precise cutting properties. The 3-ring-handle gives you a secure hold and transmits every movement reliably to the snare.

The power connection is mounted independent from the 2-ring conductor on the distal end of the handpiece. This prevents the connected power cord from moving when the snare is opened and closed, which facilitates its handling even more.

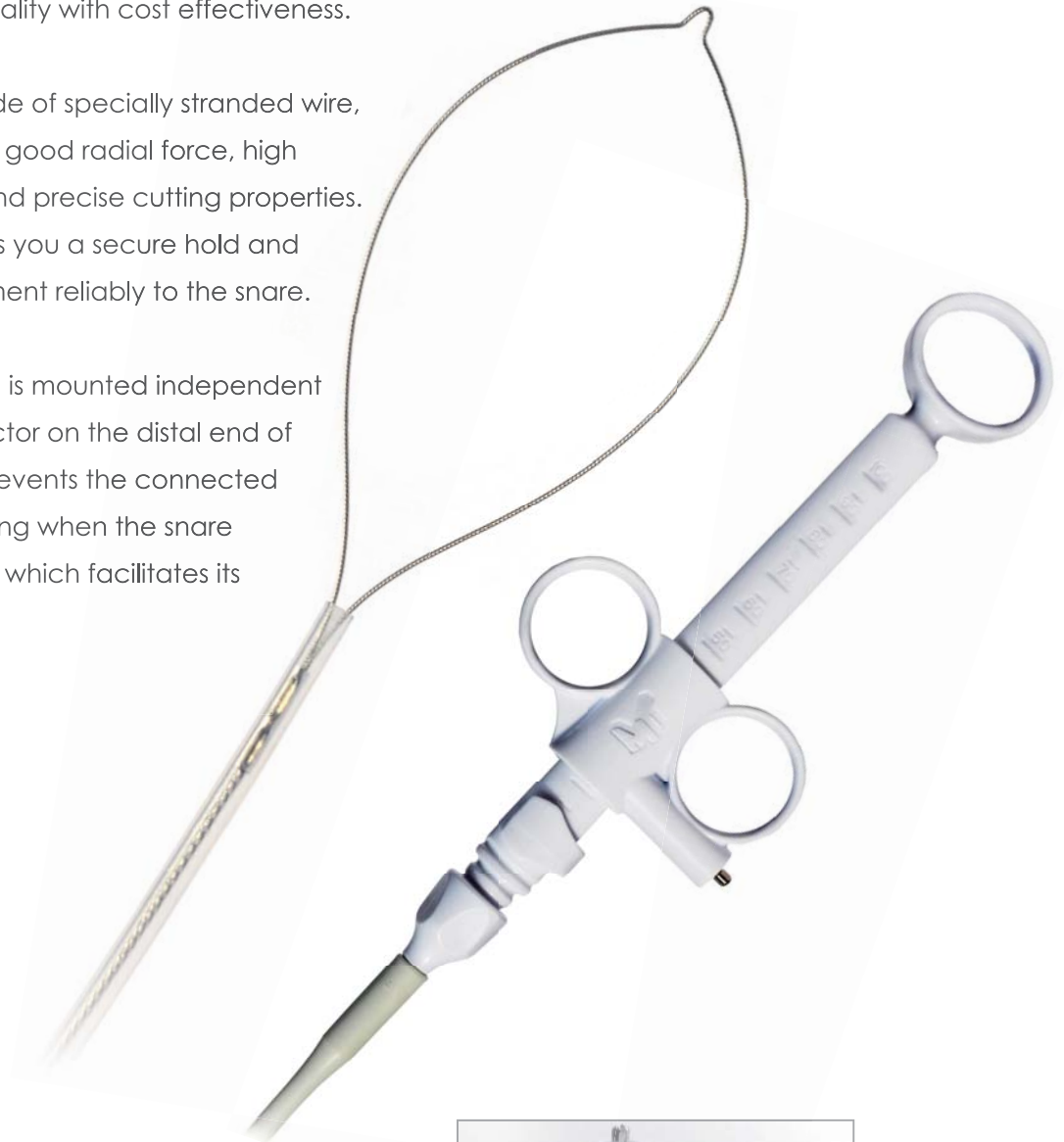


Performance Characteristics:

- very comfortable to use
- reliable and safe
- innovative power connection
- ergonomic 3-ring handle
- excellent cutting characteristics
- high strength of shape



3-ring handle with power connection

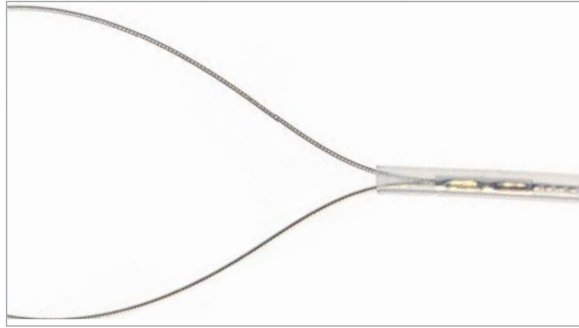

Easy and safe operation

The special stranded wire of DailySnare endows the snare with the radial force it needs to open and be controlled in any position, easily and reliably. DailySnare has proved its worth as a reliable solution in daily routine work owing to its enormous dimensional stability and good cutting properties.



The snare offers distinctly more operating comfort thanks to a simple as congenial idea: The power connection has been separated from the 2-ring conductor and can swivel around its own axis. This allows for a smooth opening of the snare even if the power cord is connected. The power cord is thus kept away from your hand.

Specifications:

REF	Technical data		Configuration	
	Ø (mm)	Length (mm)	Ø (mm)	Form
PFS01-01023230	2.3	2300	10	oval
PFS01-01523230	2.3	2300	15	oval
PFS01-02423230	2.3	2300	24	oval
PFS01-03623230	2.3	2300	36	oval