
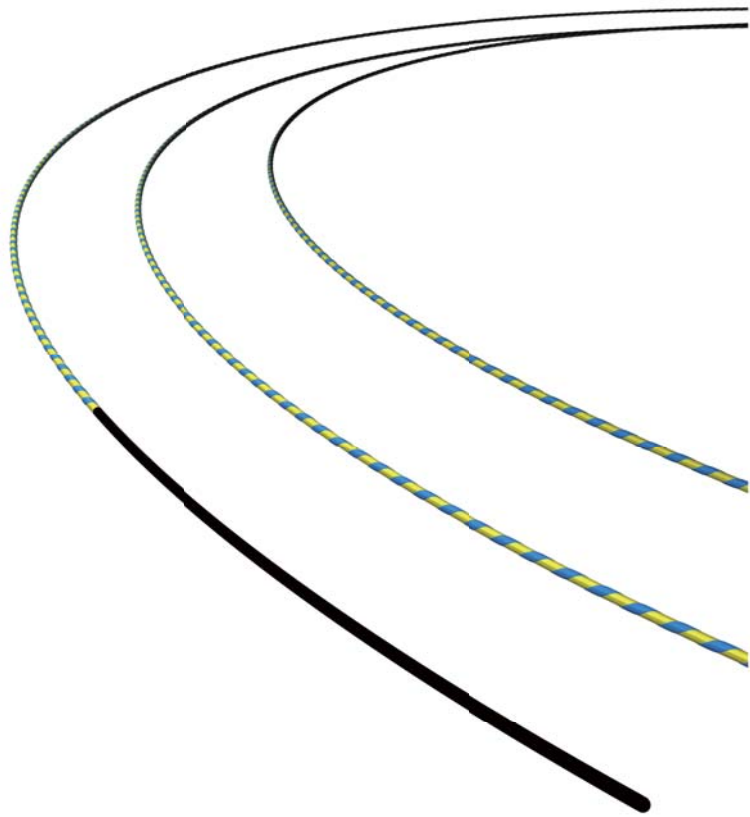


# Guide Wires

The Guide Wires program ideally complements the range of stents and thus facilitates their exact positioning and placement. All guide wires are manufactured out of bending-resistant and torsion-proof Nitinol. Owing to their hydrophilic tips, the wires safely find their way even into areas and stenoses which are hard to reach. This is supported by the wire's high rigidity and controllability.

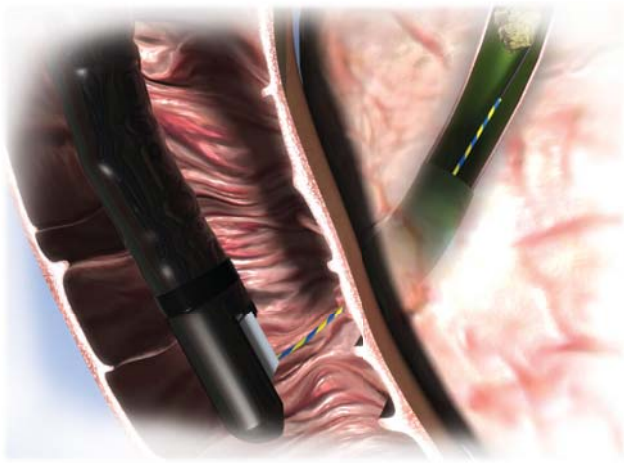

## Performance Characteristics:

- The core is made of super elastic alloy to prevent kinking.
- The atraumatic distal tip is made of soft polymeric material with hydrophilic coating to facilitate ease of cannulation.

## Specifications:

REF	Diameter (in)	Working Length (mm)	Tip Type	Color
MTN-BM-89/45-A	0.035	4500	Straight	Yellow&Blue
MTN-BM-63/45-A	0.025	4500	Straight	Yellow&Blue
MTN-BM-89/45-A-J	0.035	4500	J Style	Yellow&Blue