


## Disposable Cytology Brushes

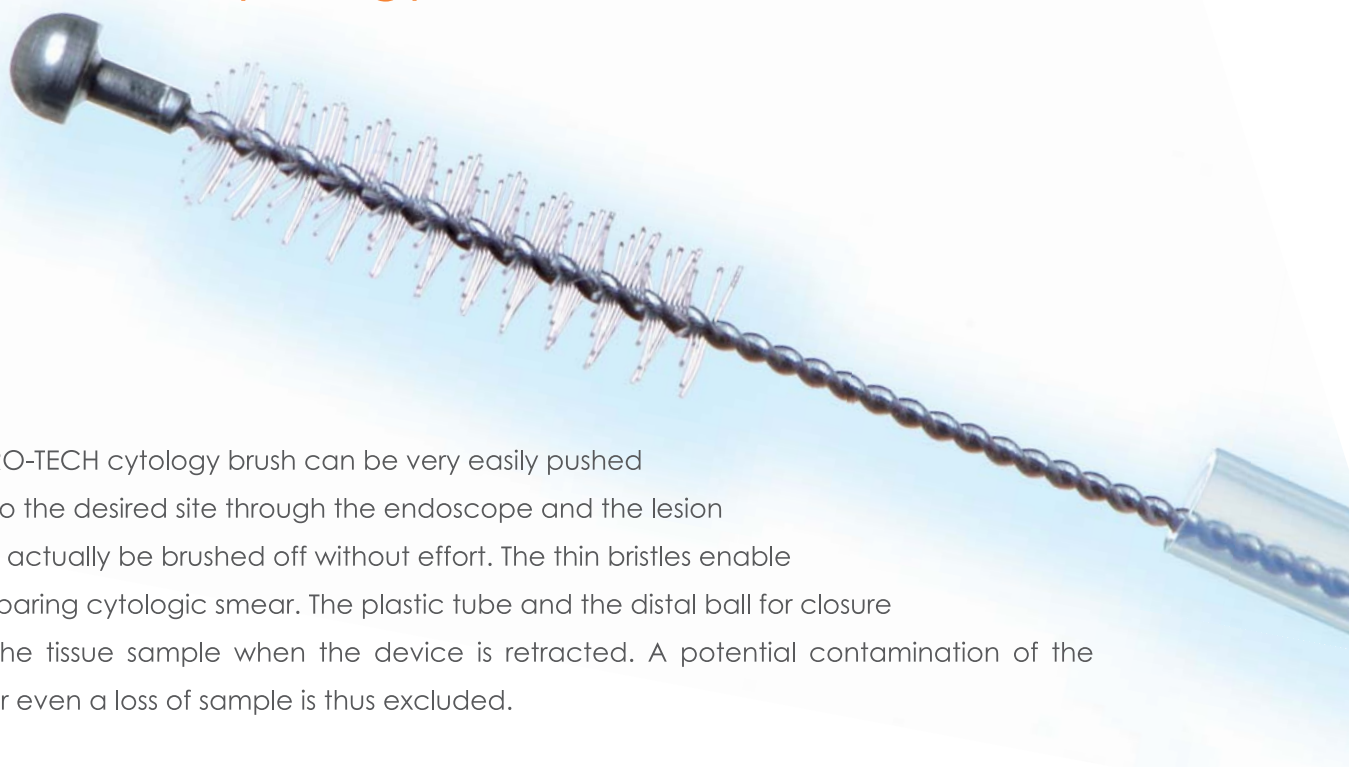

The MICRO-TECH cytology brush can be very easily pushed forward to the desired site through the endoscope and the lesion can then actually be brushed off without effort. The thin bristles enable a tissue-sparing cytologic smear. The plastic tube and the distal ball for closure protect the tissue sample when the device is retracted. A potential contamination of the sample or even a loss of sample is thus excluded.

### Performance Characteristics:

- thin brushes for an optimal collection of cells
- includes plastic tube and metal head for closure
- for endoscopes with a working channel diameter above 2.0 mm
- sterile single packaging



Ergonomically shaped 3-ring handle

### Specifications:

REF	Brush Dia. (mm)	Working Length (cm)	Catheter Dia. (mm)
CYB-24S-12C	2.4	120	1.8
CYB-30S-12C	3.0	120	1.8
CYB-30S-24C	3.0	240	2.3
CYB-40S-12C	4.0	120	1.8
CYB-40S-24C	4.0	240	2.3